




LABORATORY
EQUIPMENT

Sand baths

2021

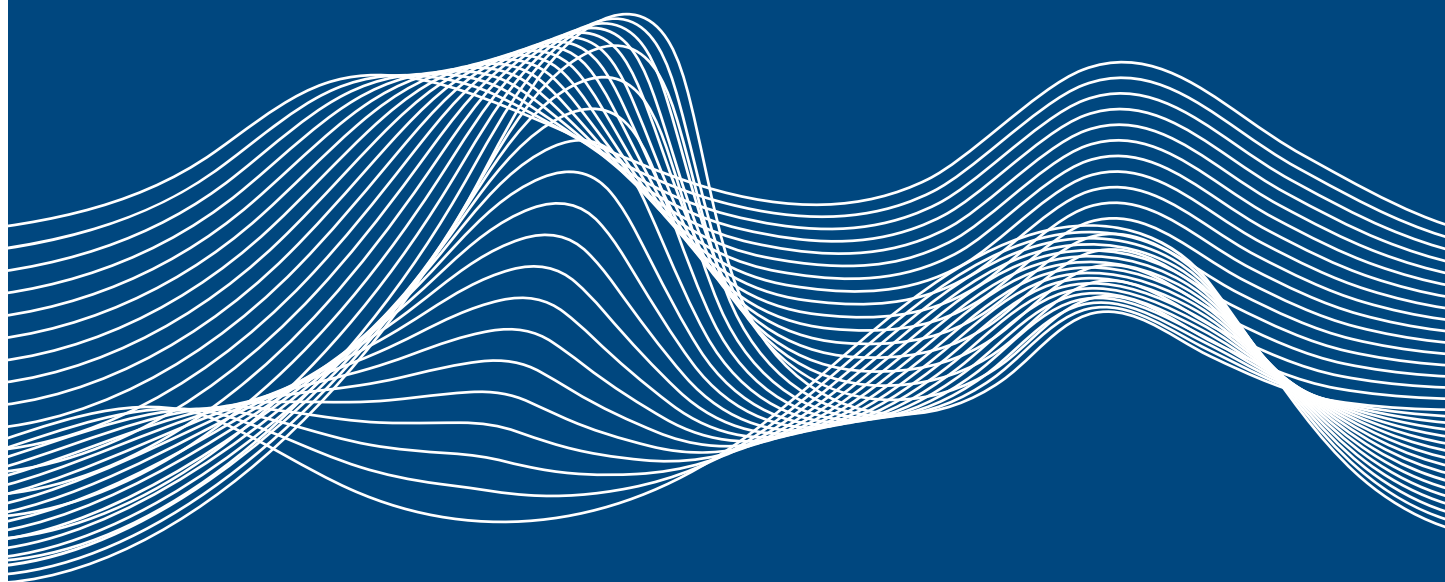

WWW.FALCINSTRUMENTS.IT

MULTIPURPOSE BATHS



Ideal for homogeneous heating of the tank



The internal tank is made of 18/8 stainless steel at high pressure



Bimetal temperature controller regulates and controls the temperature. Heating is ensured by armored heating element in stainless steel

Sand bath



S70
S130

Technical data		S70	S130
Capacity	lt	6	12
Heating power	W	1000	1400
Precision	°C	5	5
Temperature	°C	+450	+450
Sand needed	kg	6	14
External dimensions LxWxH	mm	410x400x270	610x460x270
Internal dimensions LxWxH	mm	290x230x90	500x300x90
Weight	kg	12,6	20,2
Code		614.1037.07	614.1037.13

V/HZ 220-240/50-60

SAFETY CLASS 0

PROTECTION CLASS IP 30



Suitable for tissues in paraffin, coming from microtome



The black anodized internal tank makes it easy to see the inside material. The external tank minimises heat dispersion



Temperature is controlled by a thermostat with sensor

Histology bath



BI

Technical data		BI
Capacity	lt	1,5
Heating power	W	250
Precision	°C	1,5
Temperature	°C	+30 to +90
Internal dimensions LxWxH	mm	330x95
External dimensions LxWxH	mm	200x60
Weight	Kg	1,6
Code		610.1030.00

V/HZ 220-240/50-60

SAFETY CLASS 0

PROTECTION CLASS IP 54

ACCESSORIES



Flat Lid
Code 610.1030.01



Sand in packs of 2 Kg
Code 614.1037.15



FALC INSTRUMENTS

ITALY

falc@falcinstruments.it

WWW.FALCINSTRUMENTS.IT